

Scanning Reimagined

# MEDIT Intraoral Scanner

- **MEDIT i600** - Digital Dentistry at Your FingerTips
- **MEDIT i700** - Unleash Your Clinic's Potential
- **MEDIT i700 wireless** - Magic Made Easy with a Simple Touch

## Solutions to Fuel Your Growth

Discover your daily 3D scan solutions to enhance your workflow, productivity and ensure precise results.



### Dental Clinic Solution

- Make your practice stand out
- Makes the scanning experience a comfortable one for both the dentist and patient.

### Scanning Experience


- Experience amazing digital dentistry with Medit Scan for Clinics.
- Our cutting-edge software technology makes scanning more convenient for you.

### ► Digital Dentistry with Medit Link & Medit Apps

Begin your **clinic's data management** with **Medit Link**. You can create and manage scan data matched with patient information in the electronic chart program and link them with treatment solutions of various workflows.

#### What you can do in Medit Link :

Clear aligners | Crown and bridges | Implant-supported restorations | Full and partial dentures Orthodontics | Implant surgical guides | Sleep apnea devices

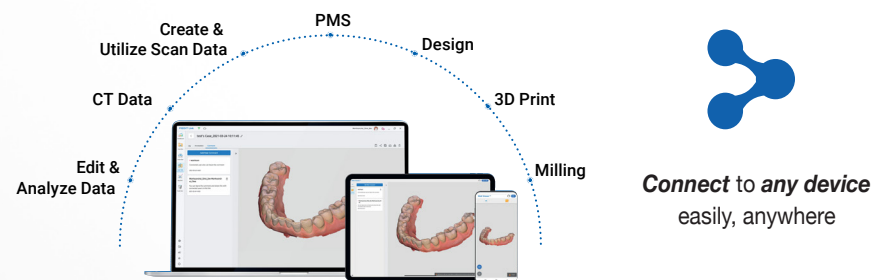

**Medit Apps** is a portfolio of software suites that puts you in the driver's seat of digital dentistry. With Medit Consultation, Medit Analysis, Medit Creator, and Medit Toolkit you are the principle architect of your workflow, adding new value to scans and images.

**Medit Consultation**

**Medit Analysis**

**Medit Creator**

**Medit Toolkit**



Category		MEDIT i600	MEDIT i700	MEDIT i700 wireless
Appearance				
Scanning Technology	Scanning Frame	Up to 35 FPS	Up to 70 FPS	Up to 70 FPS
	Imaging Technology	3D-in-motion video technology / 3D full color streaming capture		
	LED Source	LED		
	Anti-fogging Technology	Adaptive anti-fogging		
Accuracy	Full-arch 10.9 µm ± 0.98			
Handpiece	Dimensions	248 x 44 x 47.4 mm		313 x 44 x 47.7 mm
	Total Weight	245 g (incl. tip)		328 g (incl. battery & tip)
	Scan Area	15 x 13 mm		
Tip	Autoclavable	100 times	150 times	
		Autoclave 121°C 30 mins (Gravity Type) Autoclave 134°C 4 mins (Pre-Vacuum Type) Autoclave 135°C 10 mins (Gravity Type)		
	Reversible Tip	Yes	Yes	Yes
Special Mode	Remote Control Mode	No	Yes	Yes
	UV-C LED Disinfection	No	Yes	Yes
	Smart Shade Guide	Yes	Yes	Yes
	HD Mode	No	Yes	Yes
macOS Compatibility*		Yes	Yes	Yes
Connectivity	Connection Type	Wired USB 3.1 Gen 1	Wired USB 3.1 Gen 1	Wireless USB 3.1 Gen 1
	Plug & Scan	Yes	Yes	Yes (Direct connection to wireless hub)
	Detachable Cable	Yes	Yes	Yes

## MEDIT i600 Digital Dentistry at Your Fingertips

Take full control of your scanning with clicks of one button. No need for numerous cables and hubs! Simply Plug & Scan using a single cable for direct connection with a PC.

### ► Astonishing Scanning Experience

#### Light,

245g featherweight with improved grip for comfort.

#### Smooth,

Improved speed and optimal scanned images deliver a smoother scanning experience. Smooth like silk, regardless of PC specs.

#### Vivid,

Transmits vivid scan data via the new optical system used in the Medit i700.

### ► The Simplest. The Easiest



### ► Easy Does It All

Medit i600 can be **directly connected to the PC** using just **one cable**, eliminating limitations of wired devices between chairs and clinics for a wider scope of treatment using scanners.



## MEDIT i700 Unleash Your Clinic's Potential


Medit i700 makes the scanning experience a comfortable one for both the dentist and patient. With both powerful hardware and intelligent software, the full-feature Medit i700 is the key to unleash your clinic's full potential.

### ► Unleash Your Clinic's Full Potential

The Medit i700 is the key to unleash your clinic's full potential, with powerful hardware and intelligent software.

### ► Scanning experience a comfortable one

The Medit i700 makes the scanning experience a comfortable one for both dentist and patient.

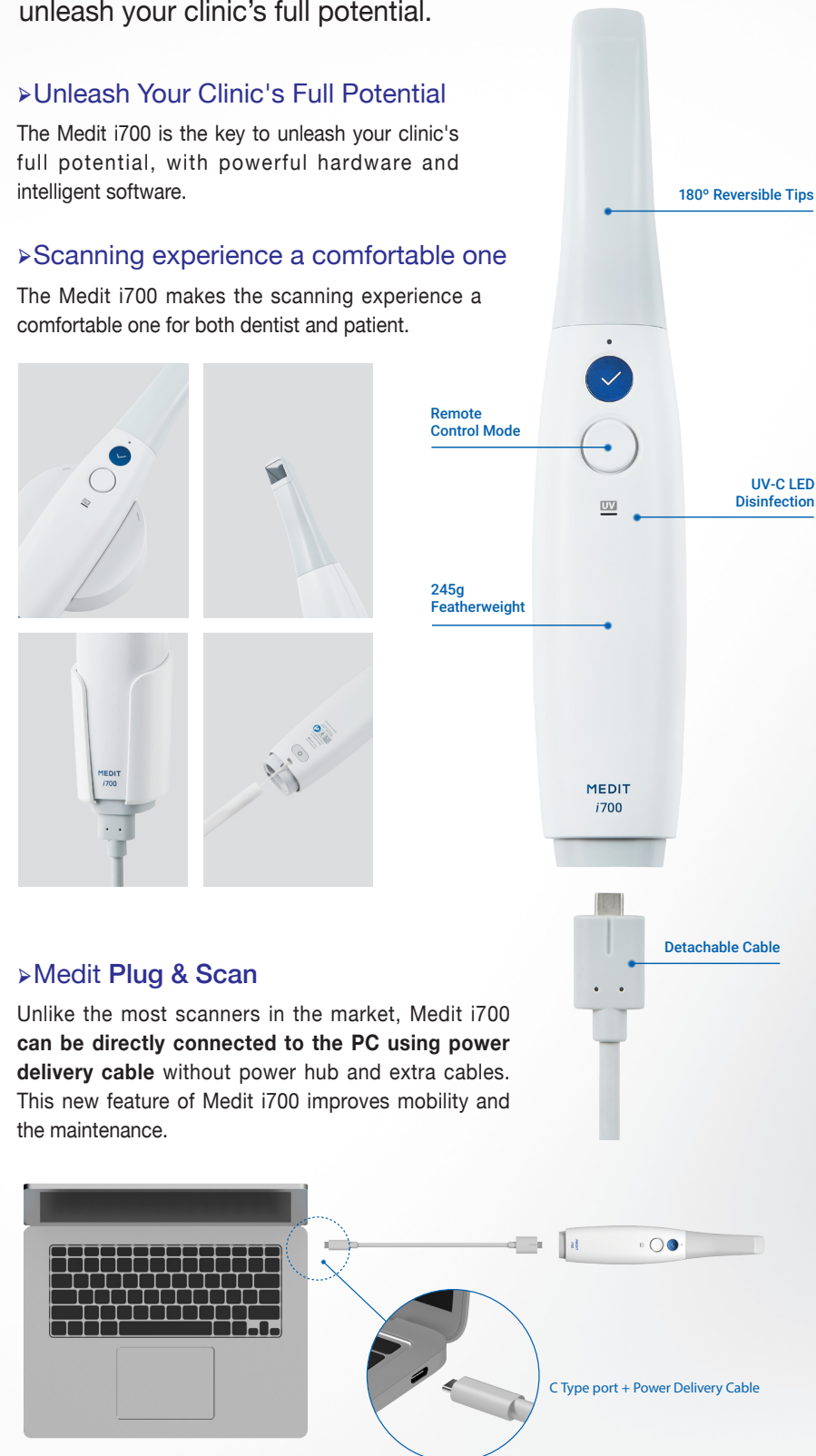

### ► Medit Plug & Scan

Unlike the most scanners in the market, Medit i700 can be **directly connected to the PC using power delivery cable** without power hub and extra cables. This new feature of Medit i700 improves mobility and the maintenance.

## MEDIT i700 Wireless Magic Made Easy with a Simple Touch

Get the power to take your practice to the next level! When the scanner's cable disappears, so do restrictions to scan movements. Strengthened impact resistance minimizes damage from accidental drops. Weight balanced design of the i700 wireless provides a stable, comfortable scanning experience

### ► No More Wires - No wire, no hassle

Free from wires, scanning is comfortable at any angle. Scan proximal areas effortlessly without constraints in movement.

### ► Beyond i700's Proven Tech

Still the same super fast, light, and accurate i700. But now **wireless**.



No More Wires



Fast, Accurate, Connected



Supercharged Yet Light



Powerful Hardware Meets Versatile Software



### ► The Best Battery Life Ever - Never Worry About Power

#### 1 1 hour of continuous scanning, up to 8 hours in standby mode.

Intelligent power management function automatically switches the device to sleep mode when not in use, standing by up to 8 hours. It takes just seconds to wake the device for scanning.



Medit i700 wireless 1 hour

Scanner A 45 min

Scanner B 30 min

#### 2 Wired Power Connection.

Once you've used up all 3 included batteries, simply connect the power cable and keep using the device.

\* Battery time can vary depending on the user environment (temperature, signal strength, etc.) and battery condition.